

Ordering information

Coronary

Distal O.D.	REF	Length	Proximal O.D.	Distal I.D.	Max GW O.D.	Marker(s)	Curve
1.7F	60172307	110cm	2.3F	0.016"	0.014"	1	Straight
	60172306	130cm				1	
	60172303	150cm				2	
	60172305	150cm				1	
	60172304	170cm				1	
	60172302	170cm				2	
1.9F	60172301	170cm	2.4F	0.017"	0.014"	3	Straight
	60192414	110cm				1	
	60192405	130cm				1	
	60192403	150cm				2	
	60192404	150cm				1	
	60192407	170cm				1	
	60192402	170cm				2	
2.6F	60192406	170cm	2.8F	0.0175"	0.014"	3	Straight Tapered
	60192401	170cm				3	
	60262801	130cm				1	
3.2F	60262802	150cm	2.3F	0.017"	0.014"	1	Dual Lumen
	60262804	170cm				1	
	60322301	135cm				2	

Peripheral

Distal O.D.	REF	Length	Proximal O.D.	Distal I.D.	Max GW O.D.	Marker(s)	Curve
2.1F	60212710	110cm	2.7F	0.021"	0.018"	1	45°
	60212711	110cm				1	Swan Neck
	60212704	110cm				1	Straight
	60212708	130cm				1	45°
	60212709	130cm				1	Swan Neck
	60212703	130cm				1	Straight
	60212706	150cm				1	45°
	60212707	150cm				1	Swan Neck
	60212705	150cm				1	Straight
	60212702	150cm				2	Straight
	60212701	150cm				3	Straight
	2.4F	60242710				110cm	2.7F
60242711		110cm	1	Swan Neck			
60242704		110cm	1	Straight			
60242708		130cm	1	45°			
60242709		130cm	1	Swan Neck			
60242703		130cm	1	Straight			
60242706		150cm	1	45°			
60242707		150cm	1	Swan Neck			
60242705		150cm	1	Straight			
60242702		150cm	2	Straight			
2.6F	60242701	150cm	3.0F	0.028"	0.025"	3	Straight
	60263001	130cm				1	Straight
2.7F	60263004	150cm	3.0F	0.028"	0.025"	1	Straight
	60273008	110cm				1	45°
	60273009	110cm				1	Swan Neck
	60273002	110cm				1	Straight
	60273006	130cm				1	45°
	60273007	130cm				1	Swan Neck
	60273001	130cm				1	Straight
	60273004	150cm				1	45°
	60273005	150cm				1	Swan Neck
	60273003	150cm				1	Straight

Please contact us for more information:





APT Medical Inc.

Address: Room 601, Building B, Tongfang Information Port, No.11 Langshan Road, Nanshan District, Shenzhen, P.R. of China

Post code: 518057 Email: international@aptmed.com

Tel: +86 755-83591013 Fax: +86 755-83480508 en.aptmed.com

PN: ENG-Microcatheter-210X285X4P-20201030

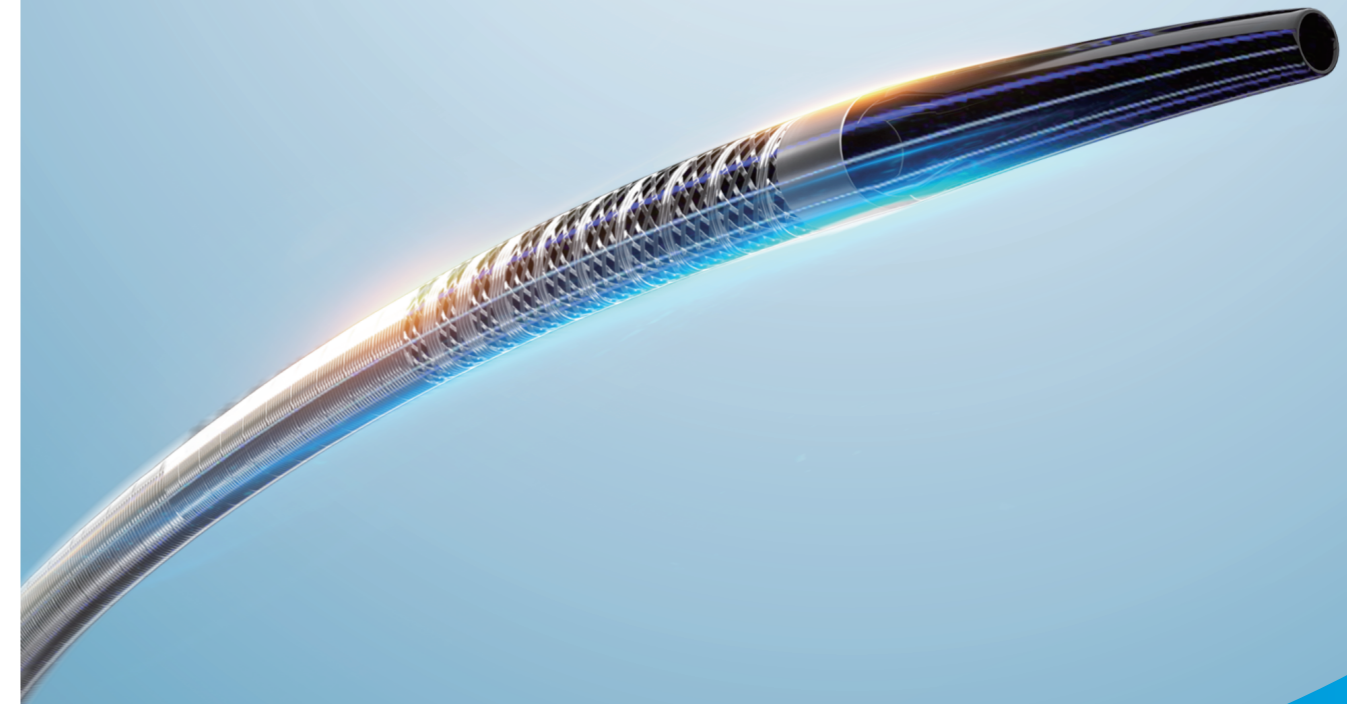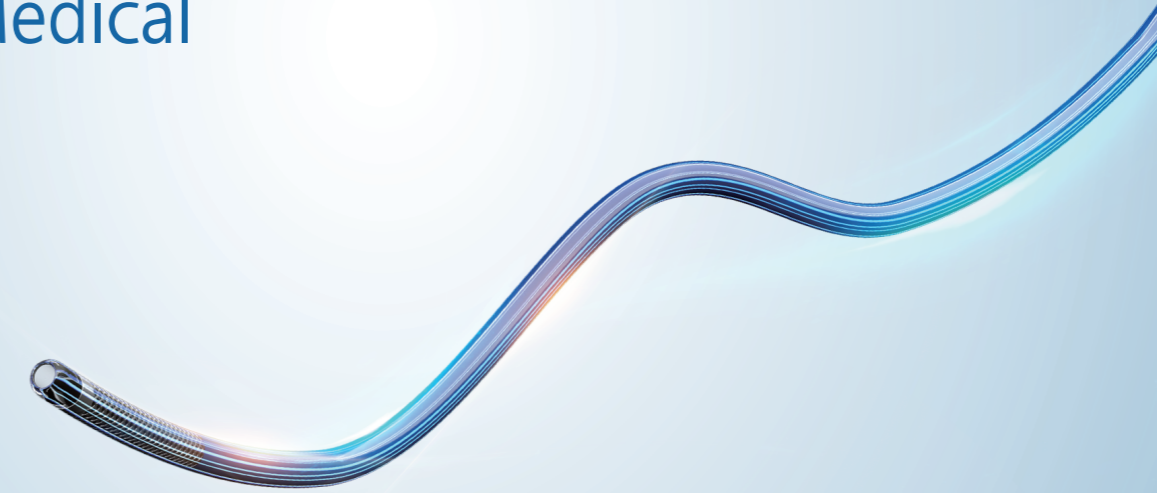

Microcatheter


Features

Optimal Support

- Stainless steel mesh braided structure provides outstanding guidewire support, superior kink resistance and device handling.

Superior resistance reduction

- PTFE inner layer for minimizing resistance during guidewire operation.
- Hydrophilic coating significantly reduces pushing resistance.

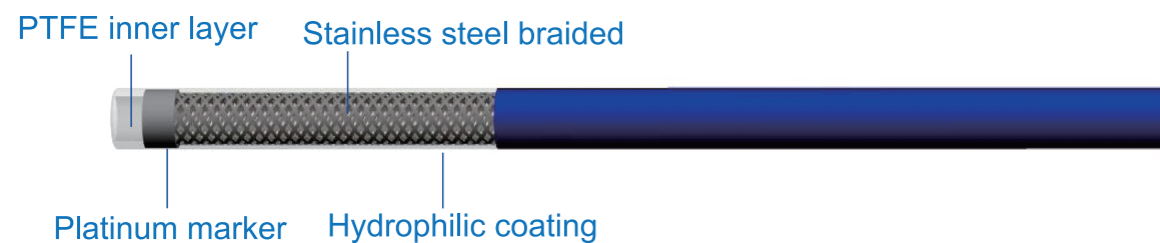


Excellent Crossability

- Multi-durometer of polymer tube offers smooth transition and excellent combination of optimal pushability, superior crossability and trackability.

Multi Radiopaque Markers

- 2 or 3 marker bands are available for easy visualization.

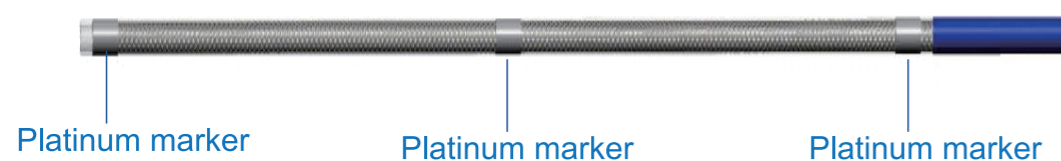

Dilatation Microcatheter



Tapered flexible tip

- For easier passing through tortuous micro channels and lesions.
- Flexible tip with tungsten ensures excellent visibility.
- Smooth transition between tip and guidewire. 0.015" tip diameter provides smooth transition for 0.014" guidewire.



Unique hybrid structure

- Reinforced coil and mesh braided hybrid structure are designed to offer excellent response, outstanding kink resistance and superior torque control.
- Provides exceptional pushability and supporting to cross challenging lesions.

Dual Lumen Microcatheter

